

PACI-TEMP THERMOMETER



MP100/MP101


PACI-TEMP BABY'S THERMOMETER

For use up to age 5 years.

Please Read Carefully Before Use:

SPECIFICATION

1. Battery: DC 1.55/1.5v (SR41/LR41)
2. Power Consumption: 0.2mW
3. Battery Life: More than 200 hours of continuous operation. (Approx. 2 years if used 10 minutes per day)
4. Display Unit: 0.1°C/ 0.1°F
5. Measuring Range: 32.0°C~42.9°C (90.0°F ~109.9 °F).
Less than 32°C/ 90°F, Lo°C / Lo°F Display
More than 42.9°C/ 109.9°F, Hi°C / Hi°F Display
6. Accuracy: ±0.1°C(±0.2°F) between 35°C~38°C (95°F~100.4°F), ±0.2°C(±0.4°F) for remainder of range at room temperature of 22°C (71.6°F)
7. Dimensions: 40mm x 43mm x 53mm
8. Weight: Approx. 14g (including battery)
9. Nipple Size: Medium

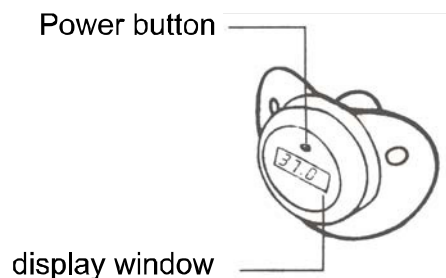

CAUTION

1. Any pacifier wears out. To PREVENT POSSIBLE CHOKING HAZARD, test the strength of the nipple by pulling on the bulb portion before each use.
2. Examine the pacifier bulb before each use for cracks, tears, tackiness or other signs of wear. Discard the worn unit immediately.
3. THE UNIT IS A MEDICAL DEVICE, do not substitute or use as a regular pacifier. Store out of reach of children.
4. Do not contact the unit with any thinner or chemical solvent.
5. Do not store the unit in direct sunlight or at high temperature.
6. Do not attempt to disassemble the unit.
7. Do not let children play with the unit. TEMPERATURE SHOULD BE TAKEN UNDER ADULT SUPERVISION.
8. When the battery is almost exhausted, the (■) symbol will appear in the display window, meaning it's time to discard the unit. Do not use as a regular pacifier.
9. ORAL TEMPERATURE READINGS ARE OFTEN (0.5°C) 1.0 °F LOWER THAN ACTUAL BODY TEMPERATURE READINGS.
10. Do not use this thermometer for measurements other than for oral use.
11. Components may be harmful if swallowed. Contact your physician immediately should this occur.
12. The Thermometer nipple comes in one size intended to fit all mouth sizes. Please gauge your child's mouth size and the nipple size prior to thermometer use. A child

using a nipple that is too large has an increase risk of choking.

HOW TO OPERATE

1. Press the on-off button to turn on the thermometer, the window will display 188.8M^oE and a beep sound will be heard. This indicates that the unit is in working condition.
2. Upon releasing the on-off button, last measured temperature will be showed for about two seconds with letter "M" under the degree sign, then 37.0°C/98.6°F, and then Lo°C/Lo°F displayed (if the room temperature is greater than 32°C/ 90°F, it will be displayed rather than Lo°C/ Lo°F). The flash °C/°F indicates the thermometer is ready to take a temperature.
3. Position the pacifier bulb into babies' mouth. Note: Drinking or eating prior to use may affect reading. Measurement should be taken one half hour after eating or drinking.
4. Once the degree sign °C/°F on the display has stopped flashing (approx. 3 minutes) and sound an alarm (Bi—Bi—Bi, if temperature exceeds 37.8°C/100°F it sounds Bi-Bi-Bi --- Bi-Bi-Bi), representing the temperature is ready to read. The reading is reflecting the oral temperature of the baby.
5. It is recommended to leave the thermometer in baby's mouth for 6~ 8 minutes to allow the nipple to warm up and get a more accurate reading.
6. The unit will automatically turn off in 10 minutes. However, to prolong battery life, it is best to press the power button to turn it off once the temperature has been taken.
7. Repeat operation procedures if interrupted during use.

CARE INSTRUCTIONS

1. Before and after each use, clean the unit with warm water and sterilize the pacifier bulb with ethyl alcohol.
2. Do not boil or wash in dishwasher.
3. Wipe clean with dry soft cloth before storage.
4. Use the original cap to store the unit. Store out of reach of children.

WARNING! DO NOT TIE PACIFIER AROUND CHILD'S NECK AS PRESENTS A STRANGULATION DANGER.